

www.ixcellsbiotech.com

10340 Camino Santa Fe, Suite C, San Diego, CA 92121 Tel: (858)412-5988 Fax: (858)368-8716 Technical Supports: <u>supports@ixcellsbiotech.com</u> Orders: <u>orders@ixcellsbiotech.com</u>

Product Information

Human Peripheral Blood CD14+ Monocytes (Positive Selection)

Catalog Number	10HU-008	Cell Number	1.0 x 10^7 cells/vial 4.0 x 10^7 cells/vial
Species	Homo sapiens	Storage Temperature	Liquid Nitrogen

Description

Monocytes are a type of white blood cells (leukocytes). Monocytes constitute 10~30% of peripheral blood mononuclear cells in the human body. They play multiple roles in immune function including replenishing resident macrophages under normal states. And in response to inflammation signals, monocytes can move quickly (approx. 8–12 hours) to sites of infection in the tissues and divide/differentiate into macrophages and dendritic cells to elicit an immune response ^[1]. CD14 is an important surface marker for monocytes.

iXCells Biotechnologies offers CD14+ monocytes isolated from normal human peripheral blood mononuclear cells (PBMCs) using positive immunomagnetic selection. > 90% of the cells are CD14+ as showed by flow cytometric analysis (Figure 1 A). And >53% of monocytes are capable to differentiate into mature dendritic cell (Fig 1B, C).

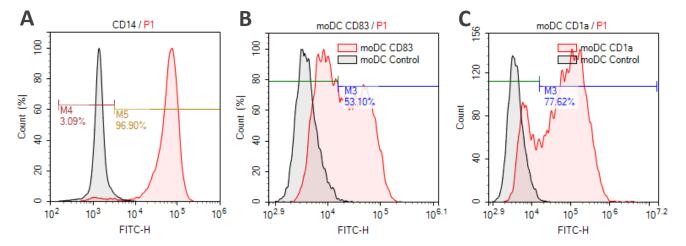


Figure 1. Flow cytometric analysis showed A. >90% of purified cells are CD14 positive. B. C. Flow cytometric analysis of monocytes dendritic cell differentiation, showed >53% of cells are CD83 positive and >77% cells are CD1a positive.

1

All Rights Reserved

Product Details

Tissue	Normal human peripheral blood
Package Size	1.0x10 ⁷ cells/vial, 4.0x10 ⁷ cells/vial
Purity	>90%
Passage Number	P0
Shipped	Cryopreserved
Storage	Liquid nitrogen
Growth Properties	Suspension
Media	Blood Cell Culture Medium (Cat# MD-0007)

Protocols

Thawing of Frozen Cells

- 1. Upon receipt of the frozen cells, it is recommended to thaw the cells and initiate the culture immediately in order to retain the highest cell viability.
- 2. To thaw the cells, put the vial in 37°C water bath with gentle agitation for 1-2 minutes. Keep the cap out of water to minimize the risk of contamination.
- 3. Pipette the cells into a 15 mL conical tube with 5 mL fresh Blood Cell Culture Medium (Cat# MD-0007).
- 4. Centrifuge at 400-450 g for 5 minutes under room temperature.
- 5. Remove the supernatant and cell is ready for downstream applications.

Safety Precaution: it is highly recommended that protective gloves and clothing should be used when handling frozen vials.

Reference

[1] Ziegler-Heitbrock, L (2007). "The CD14+ CD16+ Blood Monocytes: their Role in Infection and Inflammation, Review". *J Leukocyte Biology* **81** (3): 584–92.

Disclaimers

2

This product is intended for laboratory research purposes only. It is not intended for use in humans. While iXCells Biotechnologies uses reasonable efforts to include accurate and up-to-date information on this product sheet, we make no warranties or representations as to its accuracy. Citations from scientific literature and patents are provided for informational purposes only. iXCells Biotechnologies does not warrant that such information has been confirmed to be accurate.

This product is sent with the condition that you are responsible for its safe storage, handling, and use. iXCells Biotechnologies is not liable for any damages or injuries arising from receipt and/or use of this product. While reasonable effort is made to insure authenticity and reliability of strains on deposit, iXCells Biotechnologies is not liable for damages arising from the misidentification or misrepresentation of cultures. © iXCells Biotechnologies 2015. All rights reserved.

For Research Only

All Rights Reserved

iXCells Biotechnologies USA, Inc.